

Sinus kit



Freer elevator



Order Code : OFE

	<p>* Membrane elevator (width 5.0mm)</p> <ul style="list-style-type: none"> - To separate and elevate the sinus membrane. - To pack the bone graft material into the sinus cavity. 	<p>Length : 171mm</p>
	<p>* Periosteal elevator</p> <ul style="list-style-type: none"> - To fracture & remove the sinus bone window. - To pack the bone graft material into the sinus cavity. 	<ul style="list-style-type: none"> - To fracture and remove the sinus bone window. - To separate and elevate the sinus membrane. - To pack the bone graft material into the sinus cavity.

Membrane separator & Bone packer



Order Code : OMSC

	<p>* Cylinder shape Bone packer</p> <ul style="list-style-type: none"> - To pack the bone graft material into the sinus cavity. 	<p>Length : 167mm</p>
	<p>* Membrane separator</p> <ul style="list-style-type: none"> - To separate sinus membrane at initial stage. 	<ul style="list-style-type: none"> - To separate sinus membrane at initial stage after bone window removal. - To pack the bone graft material into the sinus cavity.



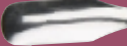

Sinus KIT
Order Code : ASLK

Sinus kit

Bone carrier with Membrane elevator





Order Code : OBGC

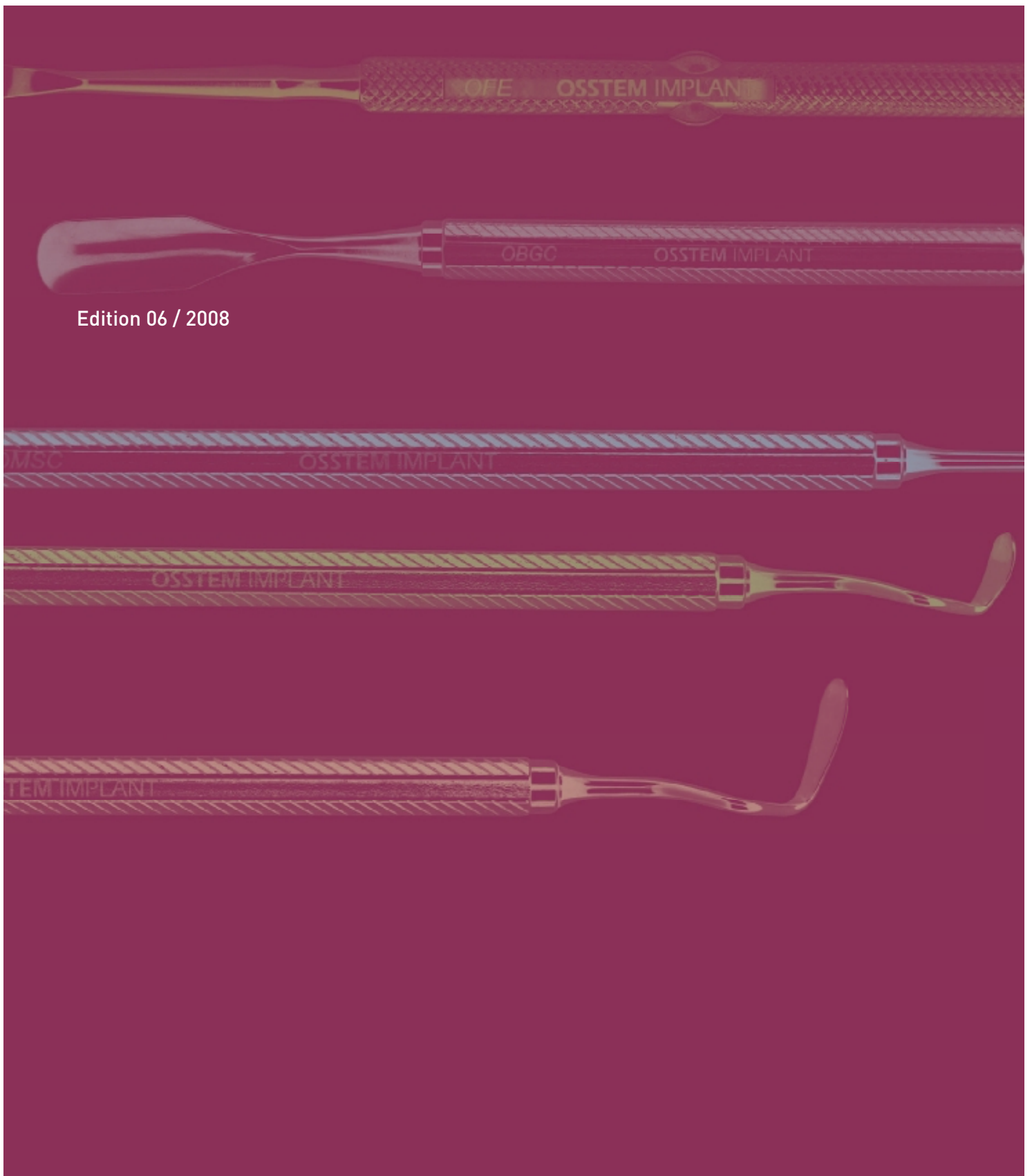
	<p>* Spatula type Bone carrier</p> <p>- To carry the bone graft material.</p>	<p>Length : 175mm</p> <p>- To carry the bone graft material on sinus lift & sinus graft technique.</p> <p>- To separate and elevate the sinus membrane through the narrow bone window.</p>
	<p>* Narrow tip Membrane elevator (width 3.0mm)</p> <p>- To separate and elevate the sinus membrane.</p>	

Sinus cures (Short & Long)



Order Code : OSCS(Short) / OSCL(Long)

	<p>* Short tip</p> <p>- Length 9mm</p>	<p>Length : Short-166mm, Long-169mm</p> <p>- To separate and elevate the sinus membrane.</p> <p>- Triple angulated tip permits easy access to any where at any directions.</p> <p>- Short or long tip can be used depend on accessibility.</p>
	<p>* Long tip</p> <p>- Length 15mm</p>	



Edition 06 / 2008

OSSTEM[®]
IMPLANT

HEAD OFFICE

8th FL. WorldMeridian Venture Center II, 426-5,
Gasam-dong, Geumcheon-gu, Seoul, Korea
tel_ 82 2 2016 7000, fax_ 82 2 2016 7001
homepage_ <http://www.osstem.com>

MANUFACTURER

#507-8, Geoje3-dong, Yeonje-gu
Busan, Korea. Zip_ 611-804
tel_ 82 51 850 2500
fax_ 82 51 861 4693